





f Cincom www.cincomhealth.com Email: service@cincomhealth.com

SHENZHEN CINCOM E-COMMERCE CO., LTD. Address: B806, building 11, gongyuandadixia oqu, Longxiangdadao ongchengjiedao,longgangqu,shenzhen, guangdong, China E-mail: service@cincomhealth.com Website: www.cincomhealth.com



TWO MODES OPTIONAL

1. Mode A: Synchronous Mode

When you press the power button, the machine will automatically work in mode A. It will provide kneading and shiatsu massage in the low air pressure intensity, you can change the air pressure intensity(L1-L2-L3)and heat(H1-H2-OF) according to your own

Feature: In this mode, the air pressure and kneading function will work simultaneously in cycle.

2. Mode B: Composite Mode

This mode provides a combination massage and will also run in the low air pressure intensity. You can change the air pressure intensity (L1-L2-L3) and heat (H1-H2-OF) according to your own preferences.

Feature: This mode will provide a different massage method. Massage functions will work in turns: air pressure massage mode-kneading massage mode-combination massage (like mode A), running in cycle.

Note: The foot shiatsu massager will be automatically shut off in 10/20/30 minutes (3 massage timings optional), you can restart it anytime if you need.

PRODUCT MAINTENANCE AND ATTENTIONS

• Do not use the device for more than 40 minutes a day.

materials which might damage the same.

• Please disconnect the cable from the power socket after use. • Ensure the device is disconnected before cleaning. • Do not clean the device with a damp cloth, brushes or hard

- Product Name: Shiatsu Foot Massager
- Model No: CM-002F Voltage: DC12V/3.0A

TROUBLESHOOTING

PROBLEM	POSSIBLE	CAUSESOLUTION
Does not work	The power supply is not on	Turn on the power supply
	Power plug is poorly connected	Check the connection
	The power plug has been damaged.	Inspect the plug for any damageConsult the professional service team for repair
	Power button has not been pressed	Press the on/off butto
A sudden interruption of the power	Power off	Check whether the power supply interface is connected well
	Massage timing is over	Press the on/off button to restart it if you need more massage

We provide 2-year warranty and long-term after-sale service. If you have any problem during using of this product, please feel free to contact us.

E-mail: service@cincomhealth.com

Note: Please write the order number together with the problems you meet in

Operating Instruction

SHIATSU FOOT MASSAGER



Thank you for choosing the CINCOM Foot massager. Please read this manual carefully and learn how to handle the device in the correct manner before using it. Pay special attention to the safety instructions and keep this manual for future reference.

Introduction —	2
Safety precautions —	2
Structure of the device —	
Operation instructions —	5
Product maintenance and attentions —	
Specification —	
Troubleshooting —	
Contact us —	

Thank you for choosing the CINCOM Foot massager. Please read this manual carefully and learn how to handle the device in the correct manner before using it. Pay special attention to the safety instructions and keep this manual for future reference.

- Nice appearance, portable and convenient.
- massage and heating function together. Kneading & shiatsu massage from toes to heels to cover all
- acupuncture points. Air pressure massage makes you more relaxed.
- 2 Massage modes: user mode and auto mode available • 3 Levels of air pressure intensities, three options for massage
- timing, two working mode available. • 2 Levels of heating temperatures adjustment: low, high
- temperature. • 10/20/30 minutes auto shut-off function.
- Powered by a DC12V adapter, safe and reliable.

SAFETY PRECAUTIONS

The product should be used in correct manners, as explained in this manual, to prevent physical and material damage.

People with any of the following conditions or those who are

- receiving medical treatments should consult the doctor prior to using the shiatsu foot masager:
- Being pregnant or at puerperium Having a fever
- The doctor has recommended you rest or you are not feeling well. Convalescing on the doctors orders.
- Suffering from a chronic illness or cancer. Having cerebral lesions (clots, lacerations, traumatism, etc.)
- Recovering from an operation to the upper part of your body, buttocks, thighs or calves.
- Suffering from a cardiac disease, heart-related problems or use an electronic medical device.
- Using a pacemaker. Injured or have a skin disease.
- Having broken or malformed bones.
- This device is not suitable for children, the disabled or invalids. • The use of the product is not recommended for people with psychic

If you feel the Air Pressure intensity too big to make you uncomfortable during use, we advise you to press the "Deflate" button for 5 seconds to relieve the discomfort.

WHILE IN USE

- Do not put the device in a place which are damp or extremely dusty, as this may result in mechanical failure.
- Do not place the device in sources of heat or in the sunlight. Place the reflexotherapy device on a flat surface. Do not use the device together with other therapeutic equipment or electrical medical apparatus and similar.
- Please avoid any liquid splashing on the equipment when using the equipment • Do not use the device on any wet part of the body, as this may result in
- mechanical failure or even an electric shock. Your body should be relaxed during treatment. • Do not use the device immediately after meals, wait for at least 1 hour. • Daily use should be limited to a maximum of 30 minutes.

Attention:

Check the voltage is suitable and never remove the cable directly.

3


WHEN THE DEVICE IS NOT IN USE


• Ensure the device is turned off after use. • Unplug the device from the mains after use and before cleaning. • Store the device in a clean and dry location.

THE DEVICE SHOULD NOT BE USED

- If it falls into water: the device should not be placed in or come into
- contact with water or other liquids. • If you feel strange or extremely uncomfortable when using the device, turn it off immediately and call your doctor.
- If there is electrical storm. • If the plug or cable is damaged or malfunctioning.

STRUCTURE OF THE DEVICE

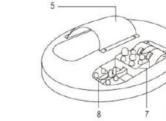

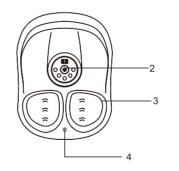

1. Power Button

3. Detachable Cloth Cover


5. Air bag 6. Front Massage Nodes 7. Arch Massage Nodes 8. Heel Massage Nodes

4. Deflate Button

2. Control Panel

Use the power adapter to connect the foot massager and the power socket, then it will make a "di" sound, and all indicators on the panel flash. The foot massager is ready for working now.

FUNCTIONS OF CONTROL PANEL

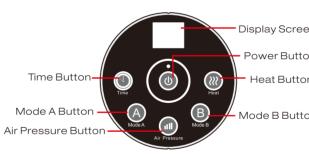

1. Press the 🗓 button to start massage(the indicator shows fixed light), the machine will run in the lowest intensity, Mode A and 20-mitnute massage timing by default.

Note: The heating function is closed by default, you can press the heat button to start it during use if you need.

2. Press 🜒 button to select massage timing: 10, 20 or 30 minutes.

3. Press button to select air pressure intensity, "L1-L2-L3" (low-medium-high) can be seen through the display screen on

4. Press button to select low or hightemperature or cancel heating: "H1-H2-OF" will be showed on the display screen on the

5.Press the 🔱 button to end the massage, the indicator will flash.