

Operating Instructions

Percussion Muscle Massager CM-029G



 Cincom

 www.cincomhealth.com

Email: service@cincomhealth.com

SHENZHEN CINCOM E-COMMERCE CO., LTD.
Address: B806, building 11, gongyuandadixiaoqu, Longxiangdadao
longchengjiedao, longgangqu, shenzhen, guangdong, China
E-mail: service@cincomhealth.com
Website: www.cincomhealth.com



RoHS
Made in China



Please read the instructions carefully before use and keep it well for further reference.

Contents

Instruction & Safety Precautions-----	01
Products Feature-----	02
Product Diagram-----	02
Operating & Charging-----	03-04
Maintenance, Cleaning & Storage-----	04
Package Included -----	04
Specifications-----	05
Contact Us-----	05

Instruction

Thanks for purchasing our product. Please read the instructions carefully before use and keep it well for further reference. If you have any problem during using of this product, please feel free to contact us through Amazon or by emailing **service@cincomhealth.com** to us. Our professional customer service team will help you!

Safety Precautions

Warnings & Cautions


- 1.Those with any of the following conditions or persons who are receiving medical treatment should consult the doctor before using the machine:
 - 1)Using pacemaker or other medical devices which are susceptible to electrical interference;
 - 2)Suffering from malignant tumors;
 - 3)Suffering from cardiac diseases;
 - 4)Having serious peripheral neuropathy dysfunction or sensory disturbance caused by diabetes;
 - 5)Being unsuitable to do the massage because of traumas on the body;
- 2.Keep it out of the reach of infants, children and people without the ability to use it independently.
- 3.Do not immerse in water Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- 4.Do not remove screws or attempt to disassemble.
- 5.Do not drop or misuse the Percussion Massage Gun.
- 6.Do not use the machine continuously for more than one hour, stop using it until you have a rest for 30 minutes.
- 7.Unplug the charger from the device after charging or before use.



8. It can only be used on soft tissues of the body, which will not cause pain or discomfort. Do not use on any hard or bony area of the head or body.
- 9.Fingers, hair, or any other part of the body should be kept away from the pole and back of the machine to avoid clipping.
- 10.Abrasions may occur during use no matter how much intensity you choose. Always check the area for treatment and stop at the first sign of pain or discomfort.

Feature

1. Powerful Brushless High-Torque Motor with Quiet Glide Technology
2. Rechargeable Lithium-Ion Battery (4-6 Hours of Use Per Charge)
3. Multi- Speed Settings (5 Speed) for Optimal Pressure
4. Cordless and Lightweight (2 Lbs)
5. ON / OFF Power Switch
6. One-button-control
7. Interchangeable Massage Heads (x6)

Product Diagram

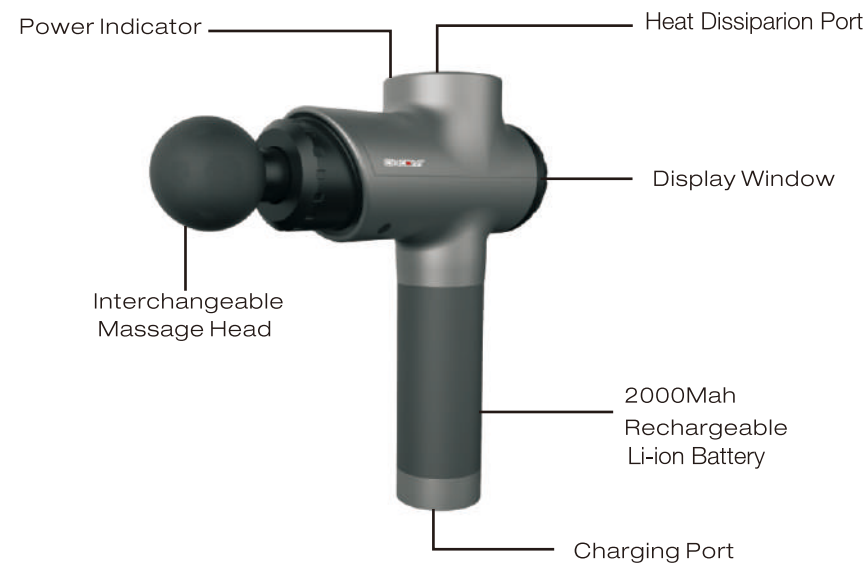

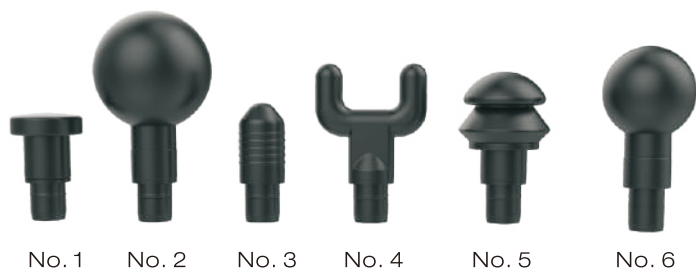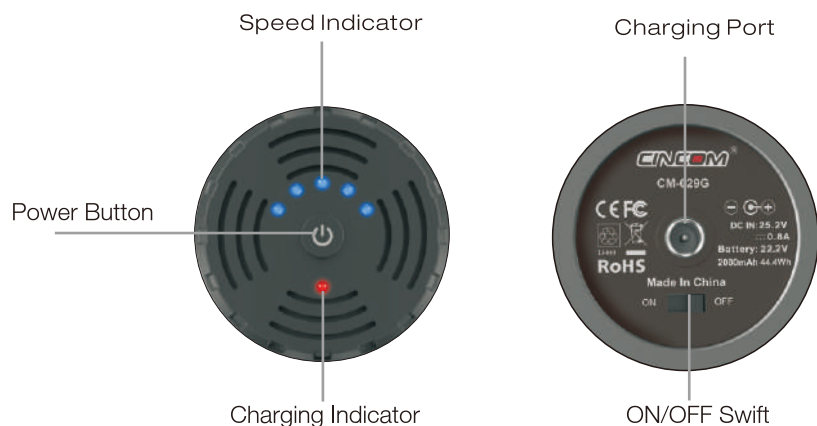

No. 1: Flat Head

Suitable for relaxing and muscle shaping.

No. 2: Big Ball Head

Suitable for back, buttocks, etc.

No. 3: Bullet Head

Suitable for meridians, palms, soles, etc.

No. 4: U-Shaped Head

Suitable for neck, spine, achilles tendon, etc

No. 5: Rubber Gas Head

Suitable for all parts of the body

No. 6: Small Ball Head

Suitable for arms, thighs, calves, etc..

Operating

1. Choose a massage head, slightly push it into the massage gun.

2. Turn on the device by moving the power switch on the bottom to "ON" position, the power indicator light is on (a blue ring light on the top of the massage gun).
3. Long press the power button(3s) to start the massage, the speed indicator light is on and it will begin with speed level 1 (the lowest intensity).
4. Adjust the speeds via the power button from Level 1-5, massage your desired body part with suitable speed, the selected speed level will be shown via the speed indicator.
5. Long press the power button(3s) to stop massage, speed indicator light will go out.
6. Turn off the device by moving the bottom switch to "OFF" position and power indicator light will go out.

Notes:

1.The device will automatically turn off after using for 15 minutes, you can press the power button to start a new session.

2.If you wanna massage with another different massage head, turn off the device, slight tug to remove the old one and slight push to put the new one in, then repeat steps 2-5.

Charging Instruction

- 1.For first charging, please make sure the machine is fully recharged, which is good for the battery.
- 2.If the charging indicator is flashing, it means that the massage gun runs out of battery. Please stop using and charge it
3. To recharge, connect the charging adapter to the device's charging port.
4. When it is recharging, the indicator of power button keeps flashing; and when it is fully recharged, the indicator is always on.
5. The average running time is 4+ hours depending on speed level and pressure applied during use.

Maintenance, Cleaning & Storage

- Use a slightly damp towel to clean device and a dry with a soft cloth.
- For storage or travel, Turn off the bottom power button and store in supplied portable box.

Package Included

- 1 x Percussion Muscle Massager
- 1 x Power Adapter/ DC 25.2V 0.8A
- 1 x Operating Instruction
- 1 x Carrying Case
- 6 x Interchangeable Massage Heads

Specifications

Model	CM-029G
Name	Percussion Muscle Massager
AC/DC Adapter	AC Input:100-240V 50/60Hz, DC Output:25.2V 0.8A
Battery Capacity	2000Mah
Working Hour	4-6H
Charging Time	2.5-3H
No-load Speed	1800-2800 Rpm
Speed Levels	5 Speeds Adjustment

Contact Us

We provide 2-year warranty and lifelong after-sale service for all products in or store, if you have any problems while using, feel free to contact us via

email: service@cincomhealth.com.

If you want to follow up our new products and latest promotional information, please visit our official website and Facebook account.